

Copyright © Vladimir Likić
PyMS version 1.0

A Python toolkit for processing of chromatography–mass spectrometry data

Contents

Introduction	1
1.1 About PyMS	1
1.2 PyMS installation	1
1.2.1 Downloading PyMS source code	2
1.2.2 Minimal PyMS installation	2
1.3 Minmax peak detector	3
1.3.1 Introduction	3
1.3.2 A brief description of the algorithm	4
Using PyMS	7
2.1 Introduction	7
2.2 PyMS examples	7
2.2.1 Example 1: Reading GC-MS data exported from Agilent ChemStation	7
2.2.2 Example 2: Exploring data attributes	8
2.2.3 Example 3: How to write an ion chromatogram or the entire intensity matrix to a file	9

Introduction

1.1 About PyMS

PyMS is a Python toolkit for processing of chromatography–mass spectrometry data. The main idea behind PyMS is to provide a framework and a set of components for rapid development and testing of methods for processing of chromatography–mass spectrometry data. An important objective of PyMS is to decouple processing methods from visualization and the concept of interactive processing. This is useful for high-throughput processing tasks and when there is a need to run calculations in the batch mode.

PyMS is highly modular and consists of several sub-packages written in Python programming language [1]. PyMS is released as open source, under the GNU Public License version 2.

There are four parts of the pyms project:

- pyms – The PyMS code
- pyms-docs – The PyMS documentation
- pyms-test – Examples of PyMS use
- pyms-data – Data used in pyms-test

Each part is a separate project on Google Code that needs to be downloaded separately.

1.2 PyMS installation

There are many ways to install PyMS depending on your computer configuration and preferences. The recommended way to install PyMS is to compile Python from sources and install PyMS within the local Python installation. This procedure is described in detail below.

PyMS has been developed on Linux, and detailed installation instructions for Linux are given below. Installation on any Unix-like system should be similar. We have not tested PyMS under Microsoft Windows.

1.2.1 Downloading PyMS source code

PyMS source code resides on Google Code servers, and can be accessed from the following URL: <http://code.google.com/p/pyms/>. Under the section "Source" one can find the instructions for downloading the source code. The same page provides the link under "This project's Subversion repository can be viewed in your web browser" which allows one to browse the source code on the server without actually downloading it.

Google Code maintains the source code by the program called 'subversion' (an open-source version control system). To download the source code one needs to use the subversion client program called 'svn'. The 'svn' client exists for all mainstream operating systems¹, for more information see <http://subversion.tigris.org/>. The book about subversion is freely available on-line at <http://svnbook.red-bean.com/>. Subversion has extensive functionality. However only the very basic functionality is needed to download PyMS source code.

If the computer is connected to the internet and the subversion client is installed, the following command will download the latest PyMS source code:

```
$ svn checkout http://pyms.googlecode.com/svn/trunk/ pyms
A    pyms/Peak
A    pyms/Peak/__init__.py
A    pyms/Peak/List
A    pyms/Peak/List/__init__.py
.....
Checked out revision 71.
$ ls -CF
pyms/
$
```

1.2.2 Minimal PyMS installation

The minimal PyMS installation involves downloading the PyMS code and installing the pycdf package. Download the pycdf package from <http://pysclint.sourceforge.net/pycdf/>, and follow the installation instructions that come with it. To test your pycdf installation:

```
$ python
Python 2.5.1 (r251:54863, Jul 10 2007, 08:23:47)
[GCC 4.1.1 20070105 (Red Hat 4.1.1-52)] on linux2
Type "help", "copyright", "credits" or "license" for more information.
>>> import pycdf
>>>
```

¹For example, on Linux CentOS 4 we have installed the RPM package 'subversion-1.3.2-1.rhel4.i386.rpm' to provide us with the subversion client 'svn'.

If the pycdf is installed properly the above import statement will produce no error messages.

The PyMS code (ie. directory pyms/) can be placed anywhere as long as the Python interpreter is made aware where it is (see Python documentation for details on how to make packages visible to Python interpreter). In example, we keep the pyms code in `/home/current/proj/PyMS/` (ie. the PyMS code is `/home/current/proj/PyMS/pyms/`). Since this is a non-standard location for Python packages, before running PyMS we set the following in Python:

```
import sys
sys.path.append("/home/current/proj/PyMS/")
```

This makes Python aware that some packages may be in `/home/current/proj/PyMS/`. Therefore to check the PyMS installation one might try the following:

```
$ python
Python 2.5.1 (r251:54863, Jul 10 2007, 08:23:47)
[GCC 4.1.1 20070105 (Red Hat 4.1.1-52)] on linux2
Type "help", "copyright", "credits" or "license" for more information.
>>> import sys
>>> sys.path.append("/home/current/proj/PyMS/")
>>> from pyms.IO.ANDI.Class import ChemStation
>>>
```

The last command has loaded the class 'ChemStation' from PyMS. Since this did not produce any error messages the PyMS was accessible (and moreover, that pycdf is properly installed, since the class 'ChemStation' uses pycdf).

For examples on PyMS use see Chapter 2.

1.3 Minmax peak detector

1.3.1 Introduction

Minmax peak detector is the simplest kind of a peak finding algoritam for TIC. It operates by finding peak maxima, and then attempting to determine peak boundaries. Cursory evidence suggests that gives results similar to the ChemStation peak detection, but this was not examined rigorously. The purpose of this algorithm is to provide an example of how 1D peak pickin can be implemented in PyMS. *At present the Minmax algorithm was not properly tested, do not use it for critical publication quality results.*

1.3.2 A brief description of the algorithm

Many peak detection algorithms are used in practice to process GC/LC-MS data, but only a few are fully documented, most notable those of open source projects MZmine [2] and XCMS [3]. MZmine detects peaks by finding local maxima of a certain width [2]. In XCMS peaks are detected by using an empirical signal-to-noise cutoff after matched filtration with a second-derivative Gaussian [3]. PyMS peak detection procedure was developed in-house, and relies on finding local maxima and local minima in the signal, followed by a subsequent refinement of peak left and right boundaries. Peak detection depends on two input parameters: window width over which a peak is expected to be a global maximum, and the scaling factor S used to calculate the intensity threshold $S\sigma$ which must be exceeded at the peak apex. The noise level σ is estimated prior to peak detection by repeatedly calculating median absolute deviation (MAD – a robust estimate of the average deviation) over randomly placed windows and taking the minimum. A detailed description of procedures for peak detection follows.

1. **Extracting local maxima.** Initially, an ordered list of local maxima in the signal with an intensity larger than a threshold is compiled. Two input parameters are specified by the user: the width of a window over which the peak is required to be a global maximum (W); and (2) the scaling factor S used to calculate intensity threshold $S\sigma$, where S is the noise level estimated previously (defaults: $W = 2$ data points, $S = 10$). User specified window is centered on each point of the signal, and the point is deemed to be a local maximum if the following is satisfied:

- (a) It is equal or greater than all of the points within the window W .
- (b) It is greater than at least one point in the half-window interval to the left, and at least one point in the half-window interval to the right.²
- (c) Any point closer to the edge of the signal than half-window is rejected.

Intensity at each local maxima is tested, and those that have the intensity below the threshold $N*S$ are rejected. Accepted local maxima are compiled into a list.

2. **Determination of peak left/right boundaries.** For each local maxima (base maximum) the stretch of the signal between itself and the next local maximum on either side is extracted. These two signal slices are searched for the first local minimum in the direction away from the base maximum point itself. The local maxima are defined in a very similar manner as the local maxima in the previous step. A point is deemed to be a local maximum if:

- (a) It is equal or smaller than all the points within the window W .
- (b) It is smaller than at least one point in the half-window interval to the left, and at least one point in the half-window interval on the right.
- (c) Any point closer to the edge of the slice than half-window is rejected. This has the effect that the boundary point cannot approach next peak's apex closer than half-window.
- (d) If no minimum point is found, set the boundary point to the point furthestest away from the base maximum, but outside to the half-window range of the adjacent peak.

²This is to reject points within intervals of uniform intensity

3. **Elimination of peak overlaps.** In spectra dense with peaks peak boundaries as found in the step (2) may overlap due to the effect of user supplied window. The list of pre-peaks is searched for overlapping peaks. In overlapping peaks the right boundary of the lower retention time peak overlaps with the left boundary of the higher retention time peak. The overlapping boundaries are resolved by finding the point of minimum intensity between the two peaks (the split point). The peak boundaries are set to one point to the left from the split point for the right boundary of the lower retention time peak, and to one point to the right from the split point for the left boundary of the higher retention time peak.
4. **Correction for long tails.** In this step peak boundaries are adjusted to remove stretches of near-uniform intensities (i.e. long tails). Each peak is divided at the apex into two halves, and each half is processed individually in the boundary-to-apex direction. A line is fitted through M points from the boundary in the least-squares sense. Prior to calculating the angle between the line and the retention time axis, the rise in intensity is normalized with the intensity at the peak apex. If this angle is below the user specified cutoff (Q) the boundary point is dropped, and the process is repeated. This adjustment is repeated until the best fit through M points from the boundary gives an angle greater than the cutoff. The parameters M and Q are user specified (defaults: $M = 3$, $Q = 1.0^\circ$).

Using PyMS

2.1 Introduction

This chapter gives examples of PyMS use. The commands used in the examples are collected into scripts and saved in the project 'pyms-test'. Each example subsection given below has a corresponding directory in 'pyms-test'. The project 'pyms-data' contains data files used in all examples. Both 'pyms-data' and 'pyms-test' are downloadable from Google Code.

The setup used for the examples below is as follows. The projects 'pyms', 'pyms-data', and 'pyms-test' were all downloaded in the directory `/home/current/proj/PyMS`. In each example the Python interpreter was made aware of the PyMS location with the following commands:

```
>>> import sys
>>> sys.path.append("/home/current/proj/PyMS/")
```

For brevity these commands will not be shown in the examples below.

Each example was run from the corresponding sub-directory of 'pyms-test': Example 1 was run in `pyms-test/example01`, and so on.

2.2 PyMS examples

2.2.1 Example 1: Reading GC-MS data exported from Agilent ChemStation

The PyMS package `pyms.IO` provides capabilities to read the raw GC-MS data stored in the ANDI-MS format. The function `IO.ANDI.ChemStation()` provides the interface to ANDI-MS data files saved from Agilent ChemStation software. The name is a reminder that this function has been tested only on the data exported from Agilent ChemStation.

ANDI-MS data format stands for Analytical Data Interchange for Mass Spectrometry, and was developed for the description of mass spectrometric data developed in 1994 by Analytical Instrument Association. ANDI-MS is essentially a recommendation, and it is up to individual vendors of mass spectrometry

processing software to implement "export to ANDI-MS" feature in their software. Furthermore, it is vendor's good will to implement ANDI-MS specifications properly. Because of these limitations one cannot be certain that one can properly read ANDI-MS files from a particular vendor without testing this first.

The file '0510_217.CDF' is a GC-MS experiment exported from Agilent ChemStation (located in 'pym-data'). This file can be loaded in the memory as follows:

```
>>> from pym.IO.ANDI.Class import ChemStation
>>> andi_file = "/home/current/proj/PyMS/pym-data/0510_217.CDF"
>>> data = ChemStation(andi_file)
-> Processing netCDF file '/home/current/proj/PyMS/pym-data/0510_217.CDF'
    [ 2784 scans, masses from 50 to 550 ]
>>>
```

The above command creates the object 'data' which is an *instance* of the class `IO.ANDI.ChemStation`. The instance 'data' has several attributes and methods associated with it.

2.2.2 Example 2: Exploring data attributes

Load the data as shown in the Example 2.

The method `get_filename()` returns the name of the file from which the data was loaded. Usage example:

```
>>> print data.get_filename()
/home/current/proj/PyMS/pym-data/0510_217.CDF
```

The method `get_ic_at_index(i)` returns an `IonChromatogram` object at index *i*. For example, to get the first ion chromatogram from the data matrix:

```
>>> ic = data.get_ic_at_index(1)
```

An `IonChromatogram` object is the one dimensional time vector containing mass intensities as a function of retention time.

An ion chromatogram is a set of intensities as a function of retention time. This can be either *m/z* channel intensities (for example, ion chromatograms at $m/z = 65$), or cumulative intensities over all measured *m/z*. In the latter case the ion chromatogram is total ion chromatogram (TIC).

The nature of an `IonChromatogram` object can be revealed by calling the method `.is_tic()`:

```
>>> print ic.is_tic()
False
```

This shows that the first ion chromatogram in the data file is for $m/z = 51$.

```
>>> tic = data.get_tic()
>>> print tic.is_tic()
True
```

This shows the method to obtain a TIC from the data. Note that both 'ic' and 'tic' are instances of the IonChromatogram class.

The method `get_ic_at_mass(mz)` returns an IonChromatogram object corresponding to given m/z . For example, to get the ion chromatogram that corresponds to $m/z = 73$:

```
>>> ic = data.get_ic_at_mass(73)
>>> print ic.is_tic()
False
```

The method `get_intensity_matrix()` returns the entire matrix of intensities:

```
>>> im = data.get_intensity_matrix()
>>> print len(im)
2784
>>> print len(im[0])
501
```

This data matrix contains 2784 time points (MS scans), and each time point corresponds to a mass spectrum of 501 m/z points.

2.2.3 Example 3: How to write an ion chromatogram or the entire intensity matrix to a file

Load the data as shown in the Example 2.

Get the TIC as shown in Example 2 and save to a file:

```
>>> ic = data.get_tic()
>>> ic.write("output/tic.dat")
```

The TIC will be written as a number of pairs (retention time, intensity), with one such pair per line:

```
$ cat tic.dat
305.666      745997
```

306.009	726566
306.352	717704
306.695	684214
307.038	701866
307.381	893306
307.724	1278099
308.067	1290984
308.410	925558
308.752	644122

[..output deleted..]

The left column is the retention time in seconds, while the right column is the corresponding intensity. This data can be conveniently loaded and plotted in matlab or octave:

```
octave:1> load tic.dat
octave:2> plot(tic(:,1)/60, tic(:,2))
```

The output is shown in Figure 2.1. In the above command the time was divided by 60 to convert the x-axis into minutes.

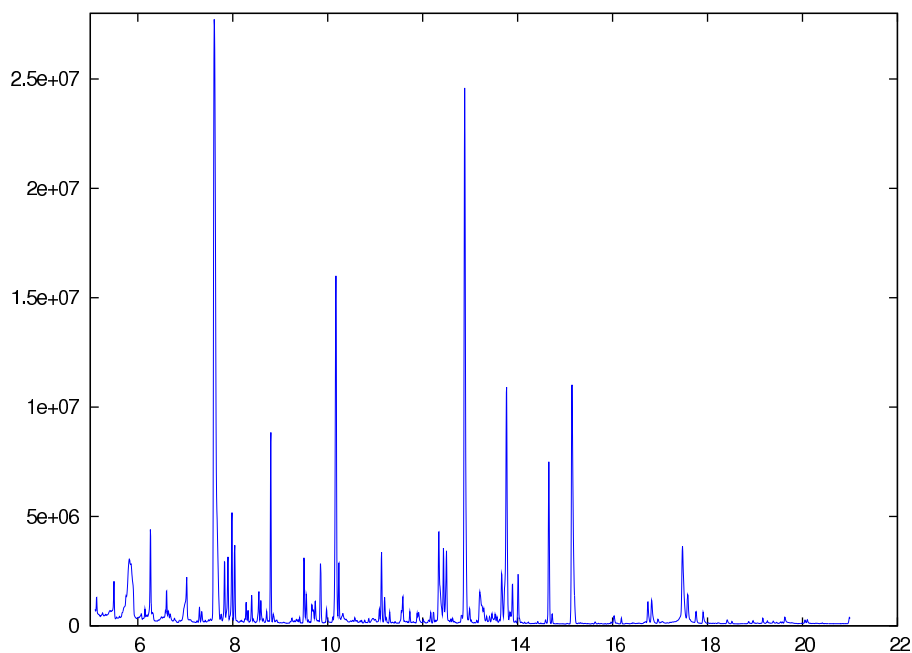


Figure 2.1: The octave plot of the file 'tic.dat'.

To get the intensity matrix and save to a file:

```
im = data.get_intensity_matrix()
```

```
save_data("output/im.dat", im)
```

The TIC can of course be calculated from the intensity matrix. For a quick check whether the intensity matrix matches the TIC one can use matlab or octave:

```
octave:1> load tic.dat
octave:2> load im.dat
octave:3> v1 = tic(:,2)';
octave:4> v2 = sum(im');
octave:6> sum(abs(v1-v2))
ans = 0
```


Bibliography

- [1] Python. <http://www.python.org>.
- [2] Katajamaa M, Miettinen J, and Oresic M. MZmine: toolbox for processing and visualization of mass spectrometry based molecular profile data. *Bioinformatics*, 22(5):634–636, 2006.
- [3] Smith CA, Want EJ, O’Maille G, Abagyan R, and Siuzdak G. XCMS: processing mass spectrometry data for metabolite profiling using nonlinear peak alignment, matching, and identification. *Anal Chem*, 78(3):779–87, 2006.

